

Topson



SoundSleep™ series CPAP and BPAP

Small in Size, Powerful Inside

Excellent User-machine Synchronization

provide better compliance.



SoundSleep machines adopt advanced flow sensor and unique flow acquisition unit to ensure the collection of small amount of flow.

Leak calculation is a key factor for good synchronization. SoundSleep machines use our patented SAF algorithm to get accurate current leak.

Finally, the breathing reversing would be finished accurately according to the volume needed.

Multilevel Respiratory Event Detection


ensures the users get the appropriate pressure.



The device can detect different flow limitation events and respond accordingly to prevent airway from collapsing in the beginning.



The device can identify different frequencies of snoring and hypopnea events accurately and adjust properly according to its severity.

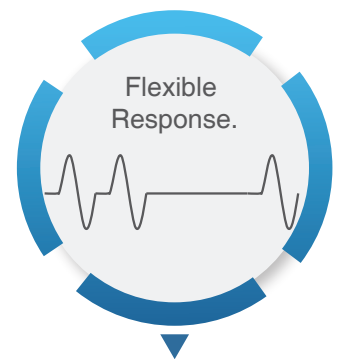

The device can detect the apnea event (OSA or CSA) and adjust the pressure flexibly.



Auto Altitude Compensation



Auto Leak Compensation



Auto on & Auto off



EPFlex



High-performance Motor



SoundSleep™ series CPAP and BPAP

User-friendly Designed Humidifier

constant temperature and humidification with multiple protection

- Open tank, Mute design
- The S-type airway ensures the smooth breathing



Easy-to-Fill
Water Chamber



Easy Clean



Metal Material



Fast Heating



Good
Humidification

Powerful Analysis Software

retrieve long-term and complete data



The Topscan software can retrieve more than 365 days data from the machines for physicians to review data and compliance. It also provides report output, print and email function.

Specifications

Parameters	CPAP	APAP	BPAP S 20/25	BPAP Auto S 20/25	S/T 20/25/30/35	S/T Pro 25/30/35	S/T TVAPS 25/30/35
CPAP	✓	✓	✓	✓	✓	✓	✓
Auto		✓	✓	✓	✓	✓	✓
S			✓	✓	✓	✓	✓
Auto S				✓		✓	✓
S/T					✓	✓	✓
T					✓	✓	✓
PC							✓
TVAPS							✓
Alarms					✓	✓	✓
Real Time Waveform						✓	✓
Pressure Range	4-20	4-20	4-20/25	4-20/25	4-20/25/30/35	4-25/30/35	4-25/30/35
Humidification	0-5	0-5	0-5	0-5	0-5	0-5	0-5
EPFlex	✓	✓					
Auto Altitude Compensation	✓	✓	✓	✓	✓	✓	✓
Auto Leak Compensation	✓	✓	✓	✓	✓	✓	✓
Noise	< 30dB						
Weight	0.95kg (without humidifier)						
Dimensions	140mm * 156mm * 94mm (without humidifier)						
Language	Multiple						

Topson Medical Technology Co., Ltd.



Address: Room 003-204, Fangda science Park, New and high-tech District, Shijiazhuang, China

Tel: 0311-67593733

Fax: 0311-67593733

Url: www.topsonmed.com